
seca 700

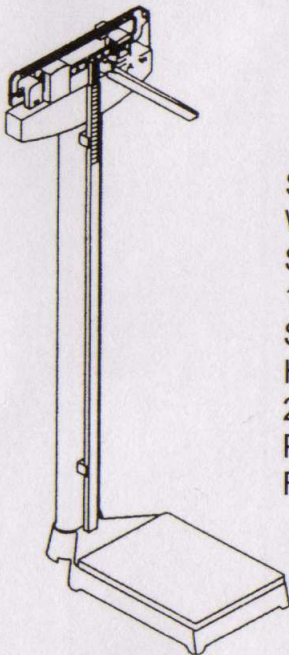SECA CORPORATION
Weighing and measuring
Systems Since 1840
1352 Charwood Road
Suite E
Hanover Maryland
21076 USA.
Ph: (410) 694 9330
Fax: (410) 694 9680

seca corporation
Model 700 112 1009
Ser. No.: 00000038972

1. Congratulations!

By purchasing the **seca** 700 personal scale you have acquired a highly accurate and sturdy piece of equipment.

For more than 150 years, **seca** has been at the service of health, and as market leader in several countries, is always setting new and higher standards with its innovative developments for weighing and measuring.

The **seca** 700 is designed for use in hospitals, doctor's offices, gymnasiums and extended care facilities.

The personal scale has a very sturdy design and will serve you a long time.

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#356

2. Safety

Before using the new scale, please take a little time to read the following safety instructions.

- Follow the instructions in the manual.
- Do not drop the scale or subject it to heavy shocks.
- Have repairs carried out by authorized persons only.

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Unpacking

- Remove all packaging

The following components are supplied:

- Weighing platform and base
- Column and head piece
- Measuring rod (inside section)
- 4 fixing screws

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
4. Assembly

- Place the column on the base with the measuring rod facing the platform.
- Screw the column to the base using 4 screws.

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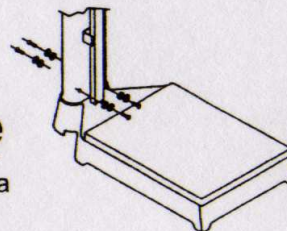
- Carefully tilt the scale to one side and hook the load lever onto the tie rod.
- Stand the scale back up again.

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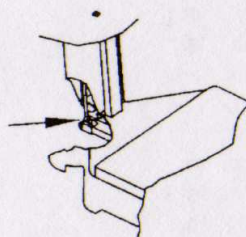
- Insert the inside section of the measuring rod into the outside section fixed to the column.
- Place the scale on a safe, level surface.

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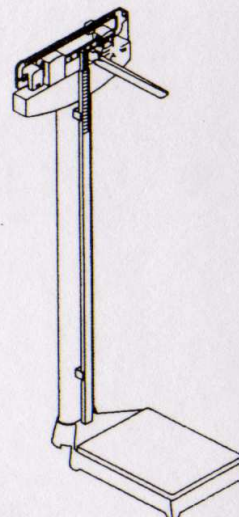
5. Setting the scale's zero point

Before the scale is used for the first time and whenever the scale is moved to another place, the scale's zero point: must be set.

First slide the two sliding weights to zero. With the scale unloaded, turn the adjuster screw until the tip of the measuring indicator is in the middle of the window.

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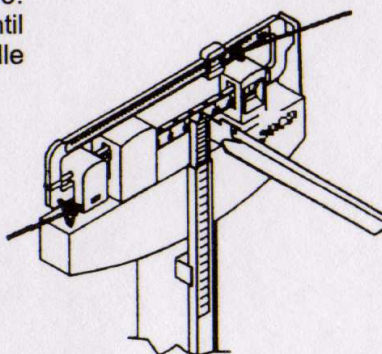
6. Operation

Correct wighing

- Step on the platform and stand still.
First of all make the coarse adjustment in steps of 20 kg. or 50 lbs.
- Lift the sliding weight on the lower scale, slide it to the right and let it latch below the roughly estimated weight.
For fine adjustment proceed as follows.
- Slide the upper sliding weight to the right until the measuring indicator settles in the middle of the window.
- Add the two values to determine the total Weight.

Fine adjustment
(in steps of - 100 g/0.25 lbs.)
Ajuste fino
(En incrementos de 100 g/0.25 lbs.)

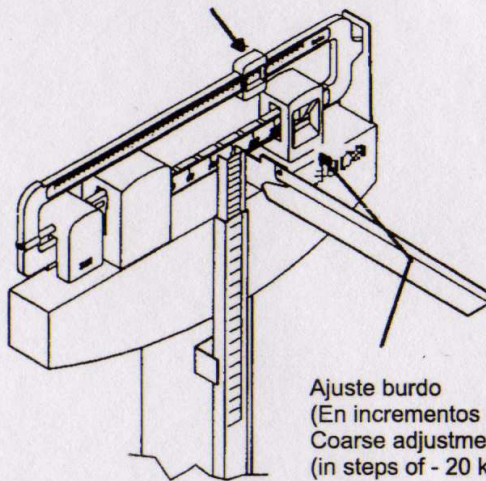Ajuste burdo
(En incrementos de 20 kg/50lbs)
Coarse adjustment
(in steps of - 20 kg/50 lbs.)

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
Measuring rod

- Fold the measuring tongue up, it will latch into the horizontal position.
- Push the measuring tongue up so that the person to be measured can easily stand below it.
- Ask the person to be measured to stand on the platform with his/her back to the measuring rod, the heels should be at the rear end of the platform; back and head must be straight.
- Now slide the measuring slide downwards until the measuring tongue rests on the head.
- For heights below 126.5 cm slide the measuring rod down until you feel resistance.
- Do not bend the measuring tongue.
- Read off the height at marking 1.

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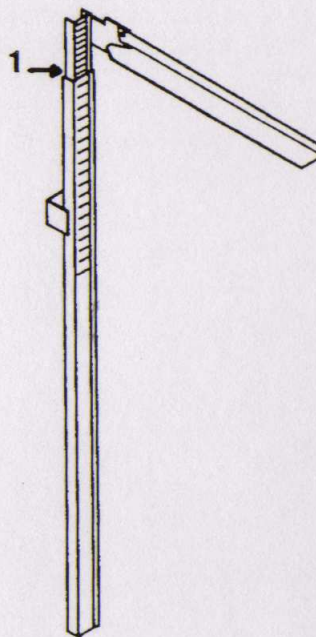7. Cleaning

Clean the outside of the scale as required with a domestic cleaning agent or a commercially available disinfectant. Follow the manufacturer's instructions.

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8. What to do if...

- the beam does not move if a person steps on the scale?
- check that the load lever is hooked in correctly.
- a weight is indicated which is obviously wrong?
- Check whether the zero point of the scale has been properly set.
- measuring errors or deviations occur when the height is determined?
- Please pay attention that the patient does not lean against the measuring rod.
- Deflection of the measuring rod will cause measuring errors.

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9. Technical data

Dimensions
Width:
Depth:
Height
Weight:
Scale:
Max. Load
Fine graduations
Measuring rod
Measuring range
Graduations

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10. Disposal

Disposal of the scale and of the measuring rod

If the scale or the measuring rod can no longer be used, your local refuse disposal service will be glad to provide information about profesional disposal in accordance with regulations.

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11. Guarantee period

We give a guarantee of 5 years from the date of purchase for defects attributable to defective materials or workmanship. All movable parts, such as batteries, cables, power packs rechargeable batteries, etc. are excluded. Defects which come under the guarantee are repaired at no charge. Other claims cannot be recognized. Transport to and from **Seca** will be at the customer's expense if the equipment is located somewhere other than at the purchaser's location. Claims under guarantee for transport damage will be honoured only if the complete original packaging was used for transportation and the scale was secured and fixed inside the packaging in accordance with its condition when originally packed. You should keep all packaging for this purpose. The guarantee is invalidated if the equipment is opened by unauthorized persons. We would ask our customers abroad to contact the sales organization of the appropriate country in the event of a claim under guarantee.

Important: the guarantee is valid only in combination with proof of purchase.

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